



STI

PCR/Syndromic panel

FluoroType[®] STI

Comprehensive STI testing from a single sample

Innovation with Integrity

Detect up to 7 sexually transmitted pathogens in one well

- Covering the most prevalent bacterial and parasitic pathogens
- Confident *Neisseria gonorrhoeae* reporting requiring dual gene detection
- Flexible workflows from low to high-throughput
- Unambiguous report providing Results at a glance

Powered by
LiquidArray[®]

IVD


Detection of sexually transmitted infections

Sexually transmitted infections (STIs) have a significant global impact, with over 1 million acquired daily and the majority being asymptomatic. Accurate and rapid diagnosis is crucial in lowering the burden these diseases present. A syndromic panel allows the comprehensive coverage of STI pathogens in a single test.

Powered by LiquidArray® | Next generation multiplexing

FluoroType® STI, powered by LiquidArray®, detects nine different targets – covering seven pathogens – in a single well. By targeting the most common bacterial and parasitic targets, FluoroType® STI offers accurate and comprehensive coverage from genital swabs and urine.

Workflow flexibility from low to high-throughput

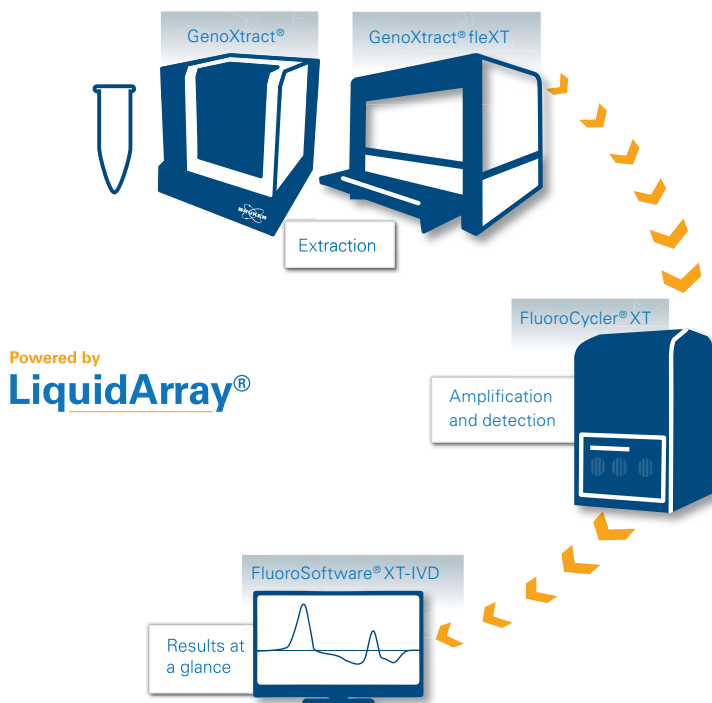

GenoXtract®: extraction 1-12 samples per run

GenoXtract® fleXT: extraction and PCR setup for up to 192 samples per day

- Maximize the number of results with each patient sample
- Minimize the number of runs per day
- Standardized reporting with Results at a glance

Replace multiple tests with one

Chlamydia trachomatis
Neisseria gonorrhoeae
Trichomonas vaginalis
Mycoplasma genitalium
Mycoplasma hominis
Ureaplasma parvum
Ureaplasma urealyticum

Confident *Neisseria gonorrhoeae* reporting

Addressing the challenge of accurate diagnosis of *N. gonorrhoeae*, FluoroType® STI detects two independent genes. The two targets are required to return a *N. gonorrhoeae* positive result, giving confidence in reporting with FluoroType® STI results.

Validated samples:

Genital swabs (endocervical, vaginal, urethral)
Urine

Results at a glance

The FluoroCycler® XT provides fully automated analysis and the integrated software generates a standardized report with results at a glance requiring minimal interpretation.

Order information



FluoroType® STI

Part No. 1877362

72 tests (GenoXtract®)

48 tests (GenoXtract® fleXT)

Please contact your local representative for availability in your country.
Not for sale in the USA, Brazil, Australia, Japan and Canada.

GenoXtract®, GenoXtract® fleXT, LiquidArray®, FluoroCycler® and FluoroSoftware® are registered trademarks of the Bruker group of companies.

GenoXtract® and FluoroCycler® XT and FluoroSoftware® XT-IVD

📍 Hain Lifescience, Hardwiesenstrasse 1, 72147 Nehren, Germany

GenoXtract® fleXT 📍 Hamilton Bonaduz AG, Via Crusch 8, CH-7402 Bonaduz, Switzerland



Hain Lifescience GmbH - A Bruker company

Nehren · Germany

Phone +49 (0) 7473-9451-0

info.mdx.de@bruker.com

hain-lifescience.de

Online information

bruker.com/microbiology

